




## **molift** PARTNER

### Technical information Molift Partner

Molift Partner - the same basic lifter for a variety of lifting tasks, gives favorable economic savings.

A full range of optional accessories quickly converts the basic lifter for use in a variety of lifting applications; sling lifting, to and from the toilet, horizontal transfers by stretcher, weighing, etc.

The Molift Partner easily copes with floor-to-bed transfers, even in intensive care units.

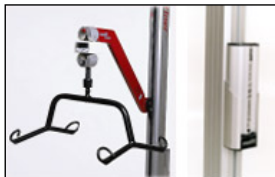

<b>Material:</b>	Lifting column in aluminum. Steel chassis and suspension
<b>Motor:</b>	12 V DC with worm gear
<b>Mechanism:</b>	Linear ball screw actuator
<b>Battery:</b>	Battery 14.4 NiCd 1.3 Ah, Molift Power pac system
<b>User load:</b>	160 kg
<b>Lifting Height:</b>	1200 mm
<b>Lifting Speed:</b>	6 double, 50 mm, polyurethane
<b>Wheels 50 mm:</b>	6 stk. doble, Ø 50 mm, polyurethanbane
<b>Wheels 115 mm:</b>	2 double, 115 mm, with brakes, rubber
<b>Legs Height:</b>	70 mm, floor to top of legs
<b>Total Width:</b>	680 mm / 800 mm with fixed legs With adjustable legs: 680 mm in inner position and 1015 mm locked in outer position (935 mm between legs)
<b>Total Length:</b>	1000 mm with lifting column mounted
<b>Total Height:</b>	1650 mm (or 1800 mm) to top of lifting column
<b>Total Weight:</b>	30 kg

### Cleaning and Maintenance

All exterior surfaces can be cleaned with regular household cleaning detergents, soapy water or hospital grade disinfectants. Abrasive compounds should be avoided. Consult your Molift User Manual for information regarding periodic maintenance.

### Key Features

- The same basic lifter for a variety of lifting tasks, gives favorable economic savings.
- The lifters are easy in use for the attendant, and comfortable and safe for the patient.
- The lifter can get under low wheelchairs and beds.
- Easy to assemble / disassemble. Easy to transport to the patients whereabouts; home care, vacation etc.
- Available with fixed / manual and electrical adjustable legs.

### Warranty:

Up to 5 years manufacturers warranty.  
1 year on sling (s), battery and hand control

### Molift Power pac

The Molift Power pac system contains of a set of batteries mounted in an aluminum profile, which are to be placed in a console mounted on the lifting column. The battery charger, a standard plug-in charger, is placed directly onto the wall.



## Technical Information

